



FILTER CHEMICAL STORAGE CABINET **HUYAIR FCS**

WWW.LABONE.VN



| UTILITIES - SAFE - EFFECTIVE |


PRODUCT OVERVIEW

FILTER CHEMICAL STORAGE CABINET - HUYAir FCS

Brand: HUYlab

Manufacturer: LABone - Vietnam

HUYlab Chemical Storage Cabinet uses activated carbon filters to help absorb the gas released during storage. The sealing system allows the cabinet to be completely sealed when closed. Each shelf has edges to prevent chemical spills and labeling. The cabinet is equipped with a central exhaust fan system, with chemical labels located. The cabinet legs can increase the height of the balance.

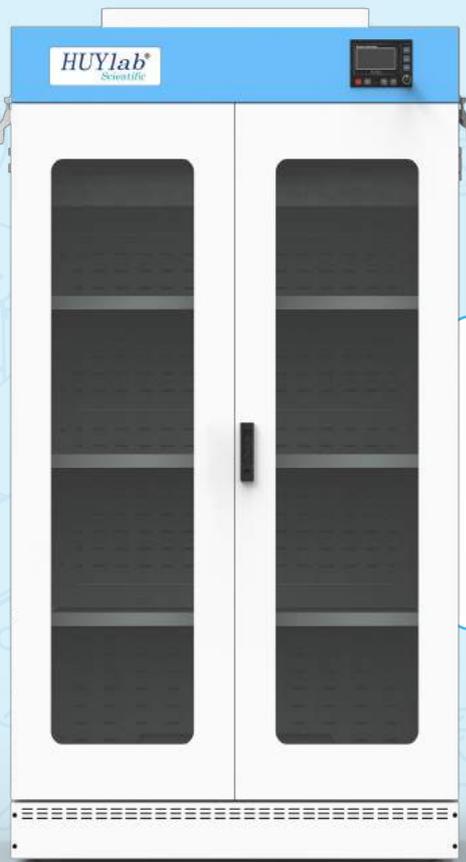


HUYlab Chemical Storage Cabinet uses activated carbon filters to help absorb the gas released during storage



PRODUCT FEATURES

- *Chemical storage cabinets with activated carbon filter help absorb a part of the gas released during storage*
- *02-layer cabinet: the entire outer structure of the cabinet is made of 1.2mm thick galvanized steel sheet with chemical resistant powder coating. Inside is 304 stainless steel material resistant to chemical corrosion*
- *The front door has a 5mm tempered glass compartment for easy viewing of the cabinet*
- *The cabinet has 4 floors for each compartment, the distance between compartments can be changed when necessary*
- *Cabinet using engine: Ebmpapst- Germany*
- *Filter size: 610x305x100 - UK*
- *The cabinet uses an LCD controller to adjust the fan speed and sensor when the cabinet is opened with strong suction airspeed to avoid the release of toxic vapors and gases. There is a pipe on the top of the refrigerator to release vapors*

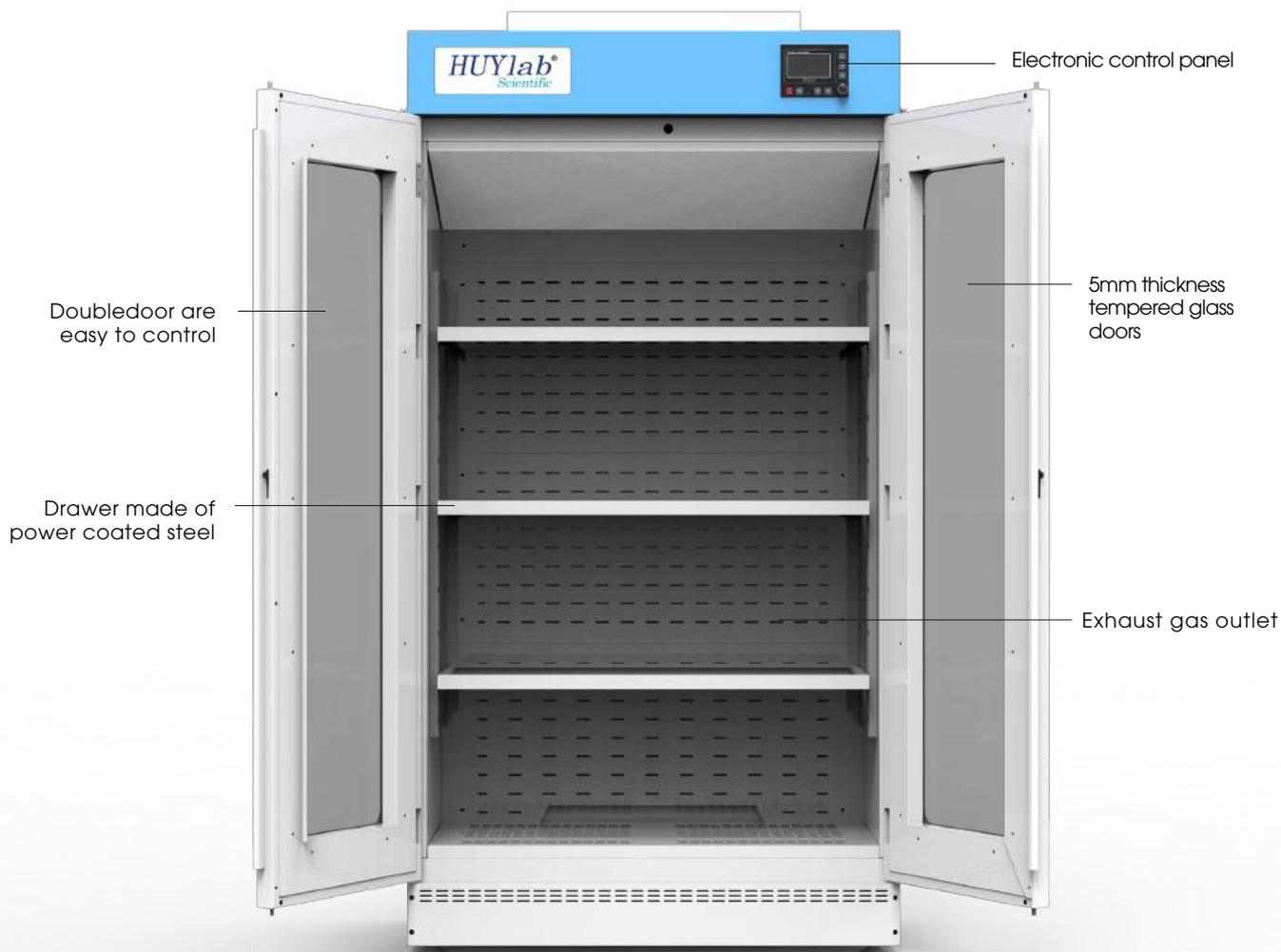

SPECIFICATION

External dimension (mm)	Internal Dimension (mm)
	850 x 550 x 2050
1200 x 600 x 2400	1150 x 550 x 2050
1500 x 600 x 2400	1450 x 550 x 2050
1800 x 600 x 2400	1750 x 550 x 2050

LABone Scientific Limited

Address: 228/13/3 Nguyen Thi Lang Street, Tan Phu Trung Communes, Cu Chi District, Ho Chi Minh City

Mail: info@huylab.com

Website: www.huylab.com